


SN450

Sterilisation instead of disinfection of instruments in the laboratory. The steriliser complies with all relevant standards and regulations for medical devices.



The indispensable safety feature for this hot air steriliser: Setpoint Wait. This means that the programme for sterilisation only starts when the set temperature has been reached. This feature can also be used with freely positionable Pt100 temperature sensors. Here the sterilisation time only begins when the set temperature has been reached at all measurement points, and reliable sterilisation is guaranteed at all times.

On this page, you can find all the essential technical data on the Memmert hot air steriliser. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at sales@memmert.com.

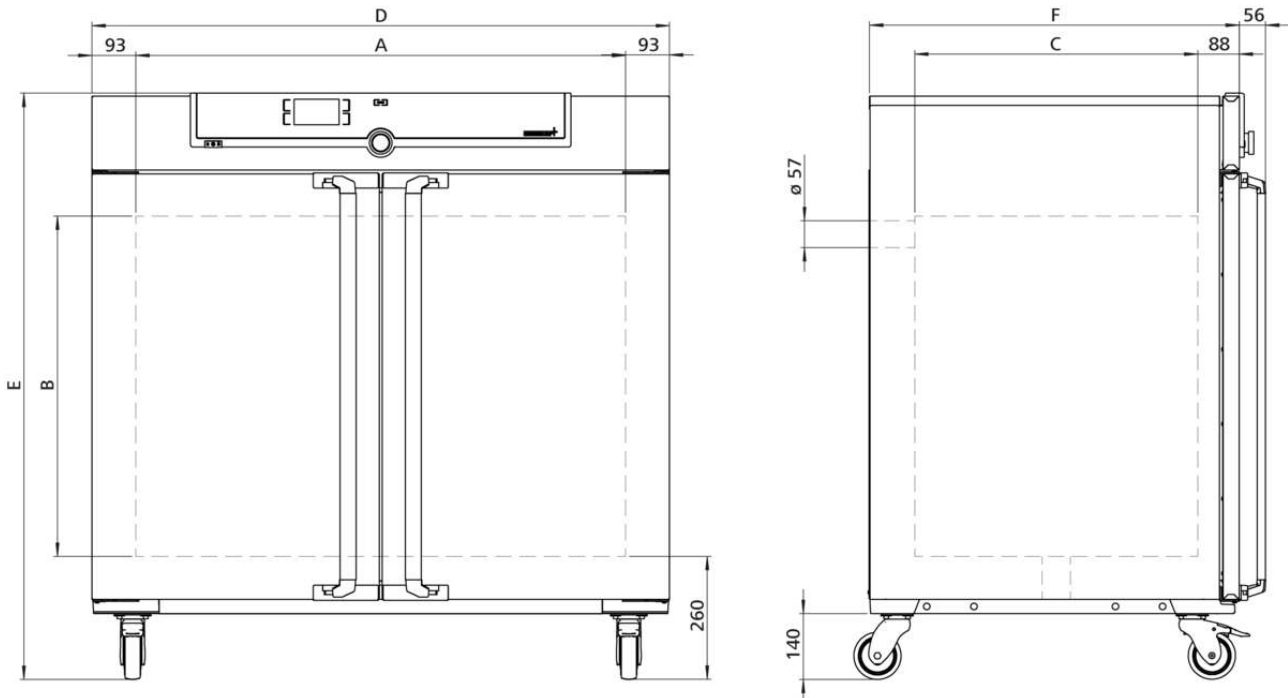

Temperature

Setting temperature range	+20 to +250 °C
----------------------------------	----------------

Setting accuracy temperature	up to 99.9 °C: 0.1 / from 100 °C: 0.5
-------------------------------------	---------------------------------------

Working temperature range	at least 5 above ambient temperature to +250 °C
----------------------------------	---

Temperature sensor	1 Pt100 sensor DIN class A in 4-wire-circuit
---------------------------	--

Control technology

ControlCOCKPIT	SingleDISPLAY. Adaptive multifunctional digital PID-microprocessor controller with high-definition TFT-colour display
-----------------------	---

Timer	Digital backwards counter with target time setting, adjustable from 1 minute to 99 days
--------------	---

Function SetpointWAIT	the process time does not start until the set temperature is reached
------------------------------	--

Calibration	three freely selectable temperature values
--------------------	--

adjustable parameters	temperature (Celsius or Fahrenheit), air flap position, programme time, time zones, summertime/wintertime
------------------------------	---

Ventilation

Convection	natural convection
-------------------	--------------------

Fresh air	Admixture of pre-heated fresh air by electronically adjustable air flap
------------------	---

Vent	vent connection with restrictor flap
-------------	--------------------------------------

Communication

Documentation	programme stored in case of power failure
----------------------	---

Programming	AtmoCONTROL software for reading out, managing and organising the data logger via Ethernet interface (temporary trial version can be downloaded). USB stick with AtmoCONTROL software available as accessory (on demand).
--------------------	---

Safety

Temperature control	adjustable electronic overtemperature monitor and mechanical temperature limiter TB, protection class 1 according to DIN 12880 to switch off the heating approx. 20°C above nominal temperature
----------------------------	---

Autodiagnostic system	for fault analysis
------------------------------	--------------------

Standard equipment

Works calibration certificate	Calibration at +160°C
--------------------------------------	-----------------------

Door	fully insulated stainless steel doors with 2-point locking (compression door lock), lockable
-------------	--

Internals	2 stainless steel grid(s), electropolished
------------------	--

Stainless steel interior

Interior	easy-to-clean interior,made of stainless steel, reinforced by deep drawn ribbing with integrated and protected large-area heating on four sides
Volume	449 l
Dimensions	w _(A) x h _(B) x d _(C) : 1040 x 720 x 600 mm
Max. number of internals	8
Max. loading of chamber	300 kg
Max. loading per internal	30 kg

Textured stainless steel casing

Dimensions	w _(D) x h _(E) x d _(F) : 1224 x 1247 x 784 mm (d +56mm door handle)
Installation	on lockable castors
Housing	rear zinc-plated steel

Electrical data

Voltage	400 V and 3x 230 V w/o neutral, 50/60 Hz approx. 5800 W
Electrical load	

Ambient conditions

Set Up	The distance between the wall and the rear of the appliance must be at least 15 cm. The clearance from the ceiling must not be less than 20 cm and the side clearance from walls or nearby appliances must not be less than 5 cm.
Altitude of installation	max. 2,000 m above sea level
Ambient temperature	+5 °C to +40 °C
Humidity rh	max. 80 %, non-condensing
Overvoltage category	II
Pollution degree	2

Packing/shipping data

Transport information	The appliances must be transported upright
Customs tariff number	8419 8998
Country of origin	Federal Republic of Germany
WEEE-Reg.-No.	DE 66812464
Dimensions approx incl. carton	w x h x d: 1330 x 1440 x 1050 mm
Net weight	approx. 161 kg
Gross weight carton	approx. 227 kg

Standard units are safety-approved and bear the test marks

